



Peptide Identification (Repetition)


Shotgun proteomics

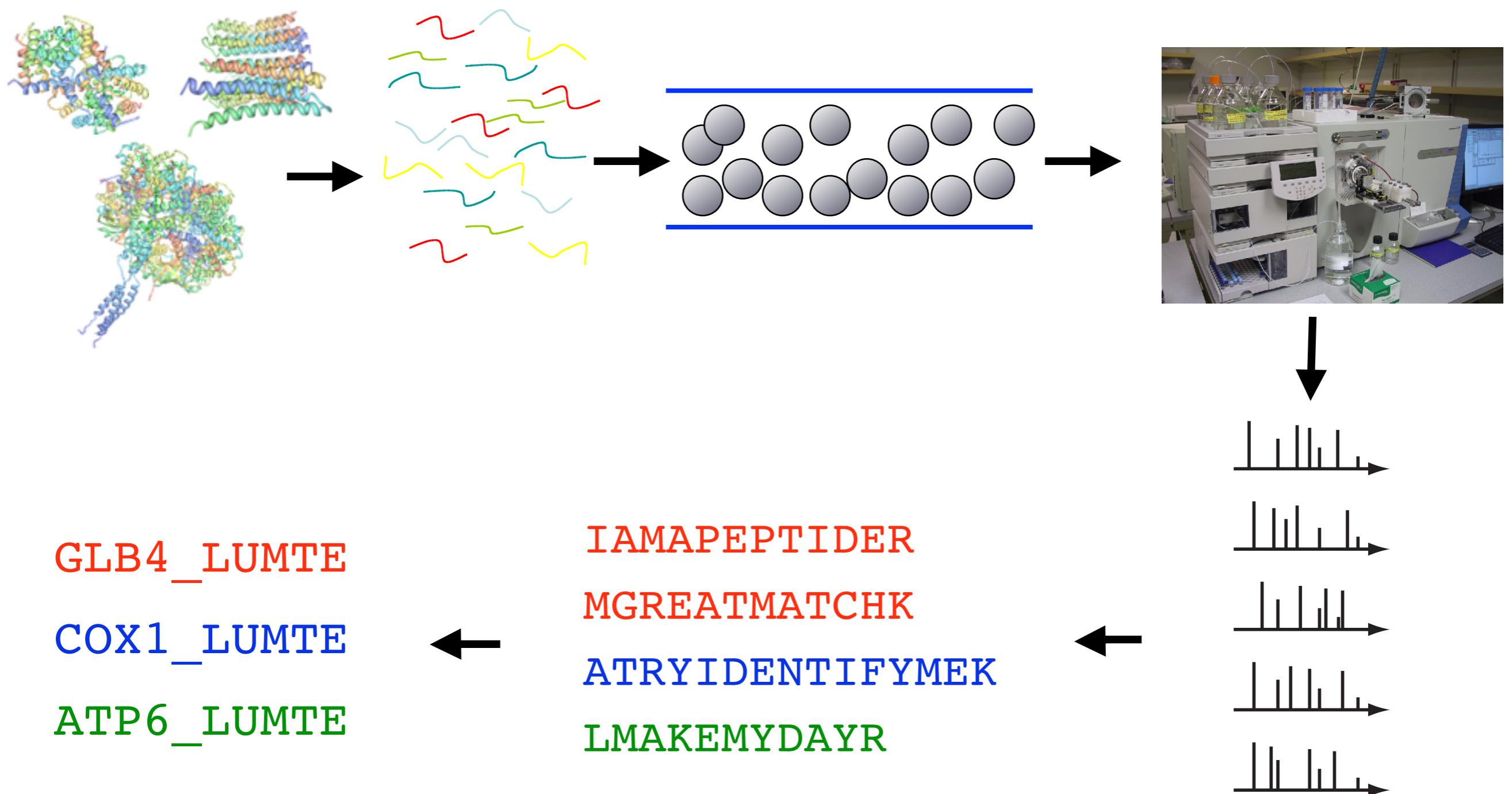


Peptide spectra

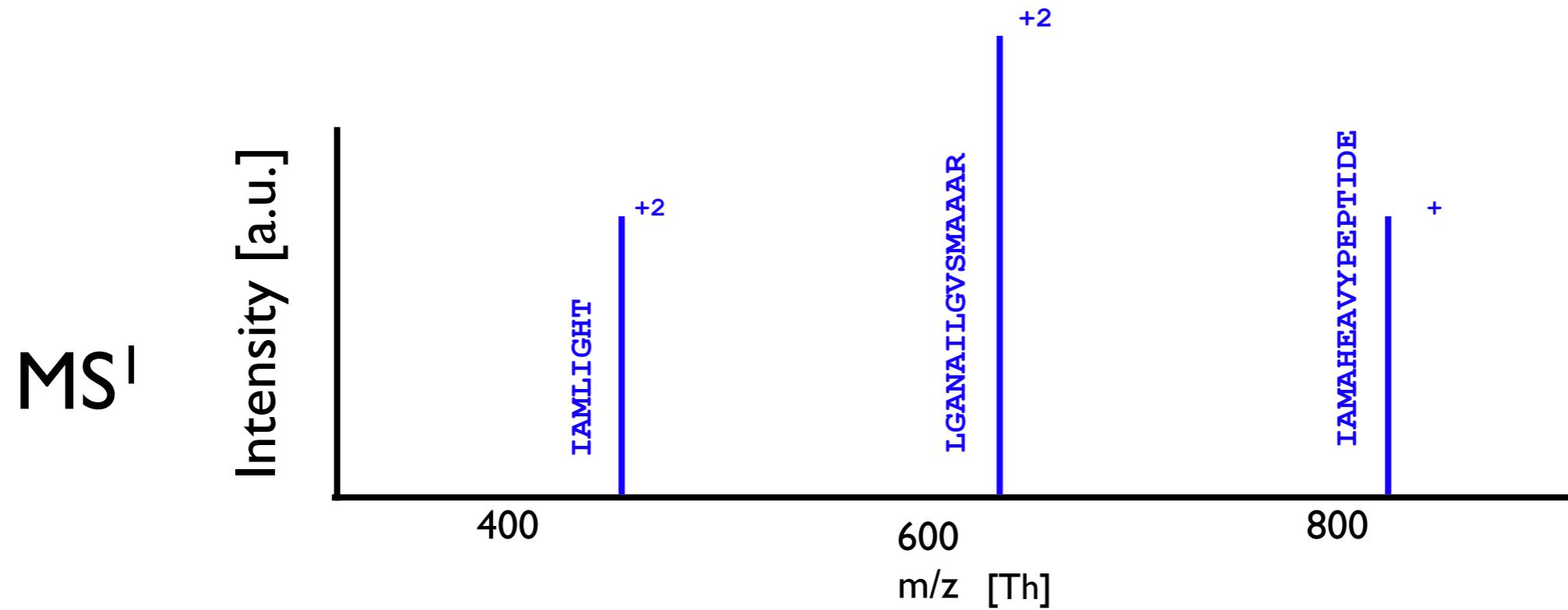



Peptide spectra

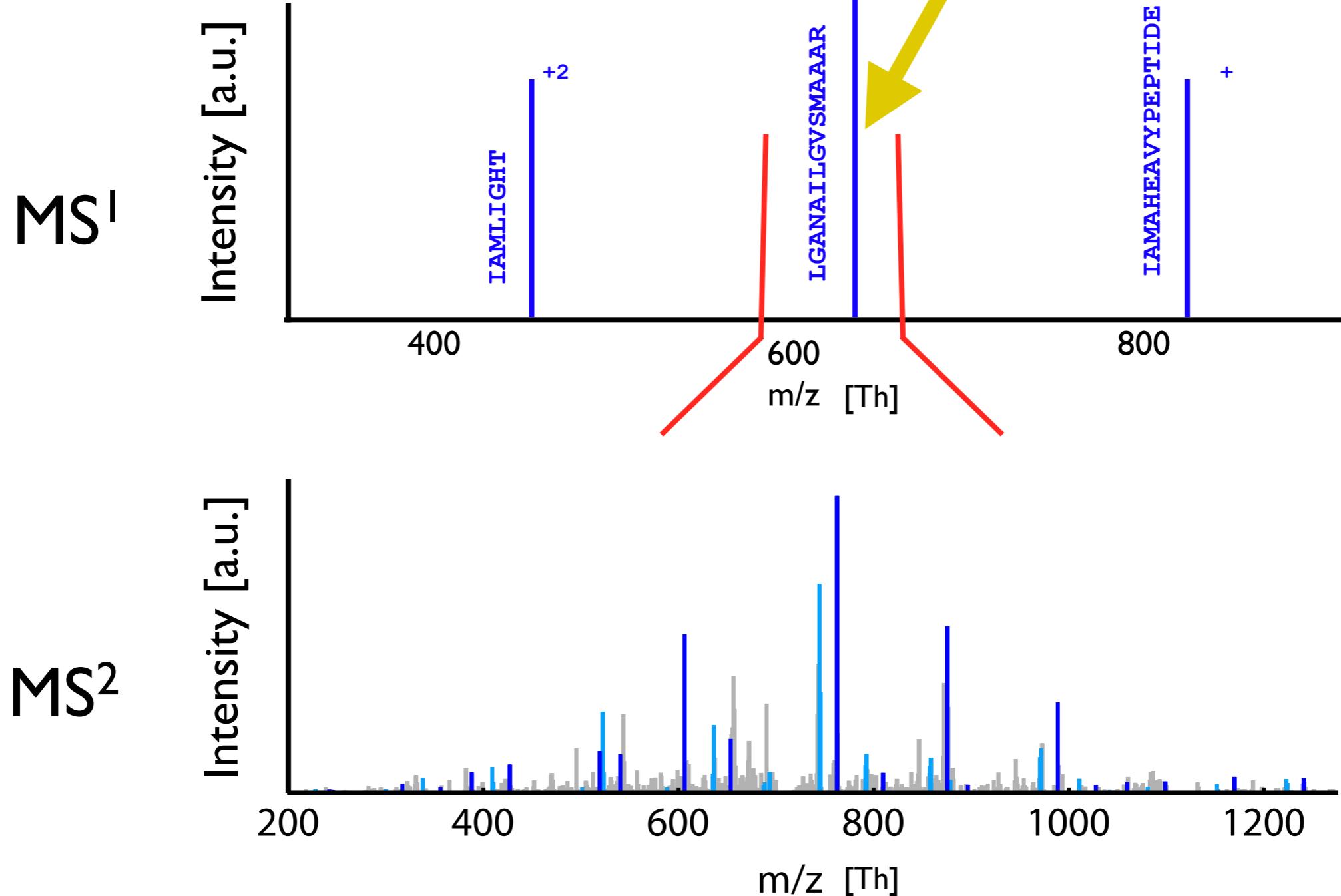


Chromatograms and Fragmentation spectra

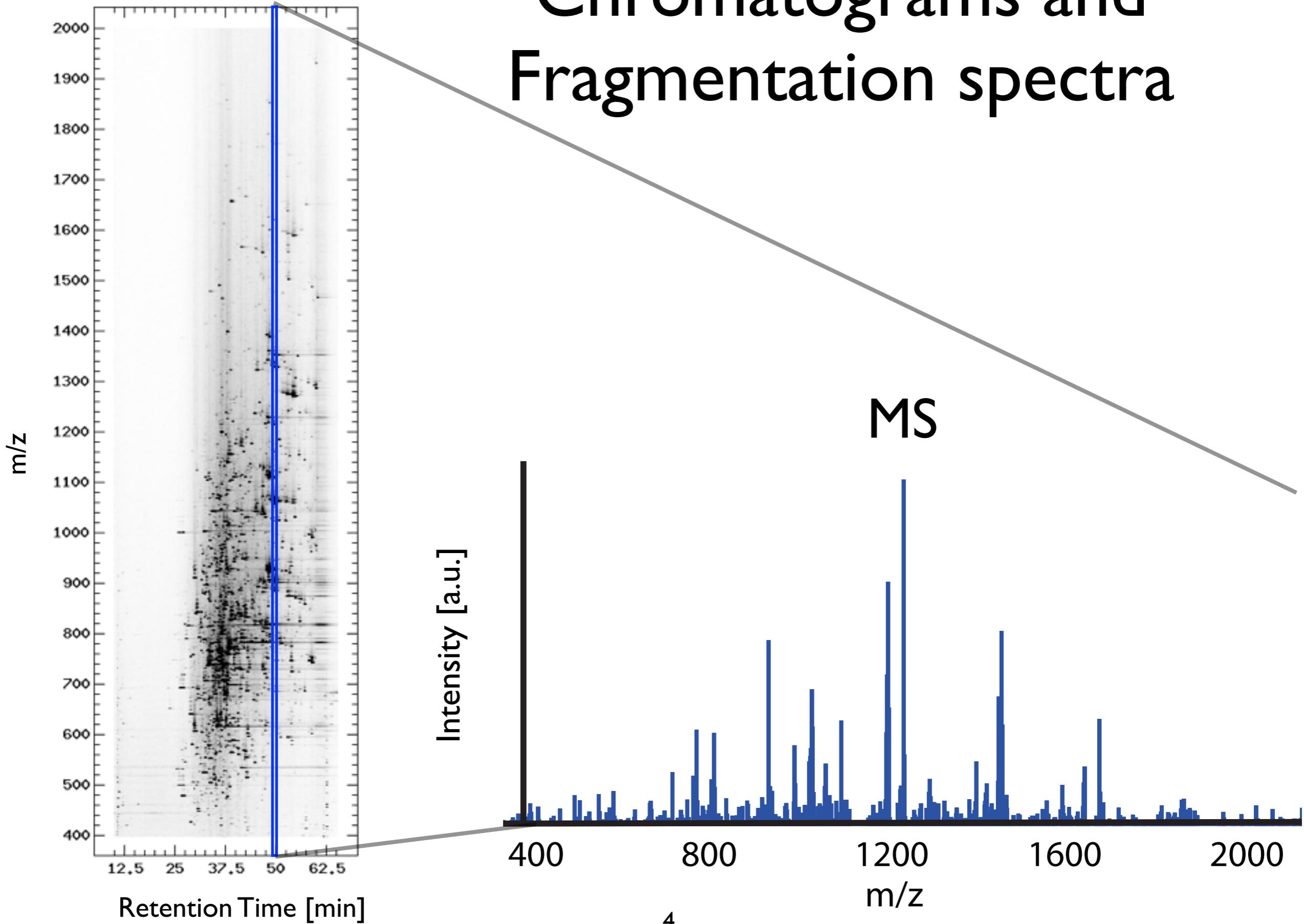



Chromatograms and Fragmentation spectra

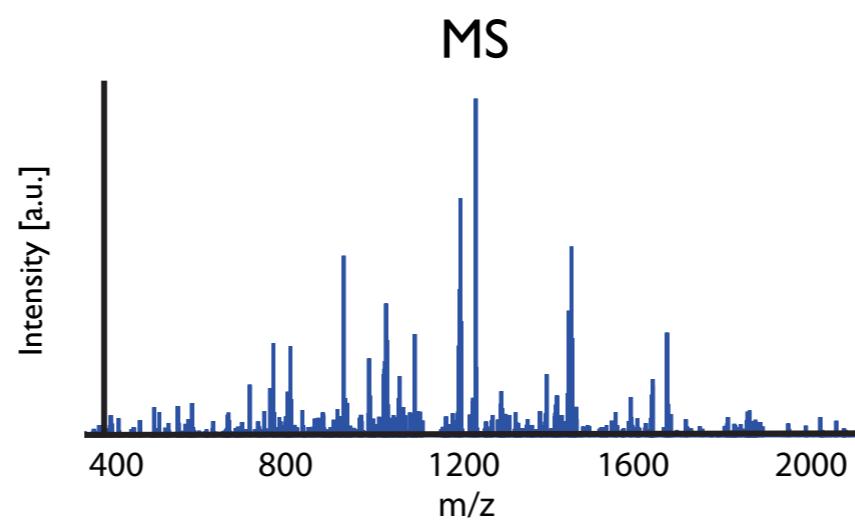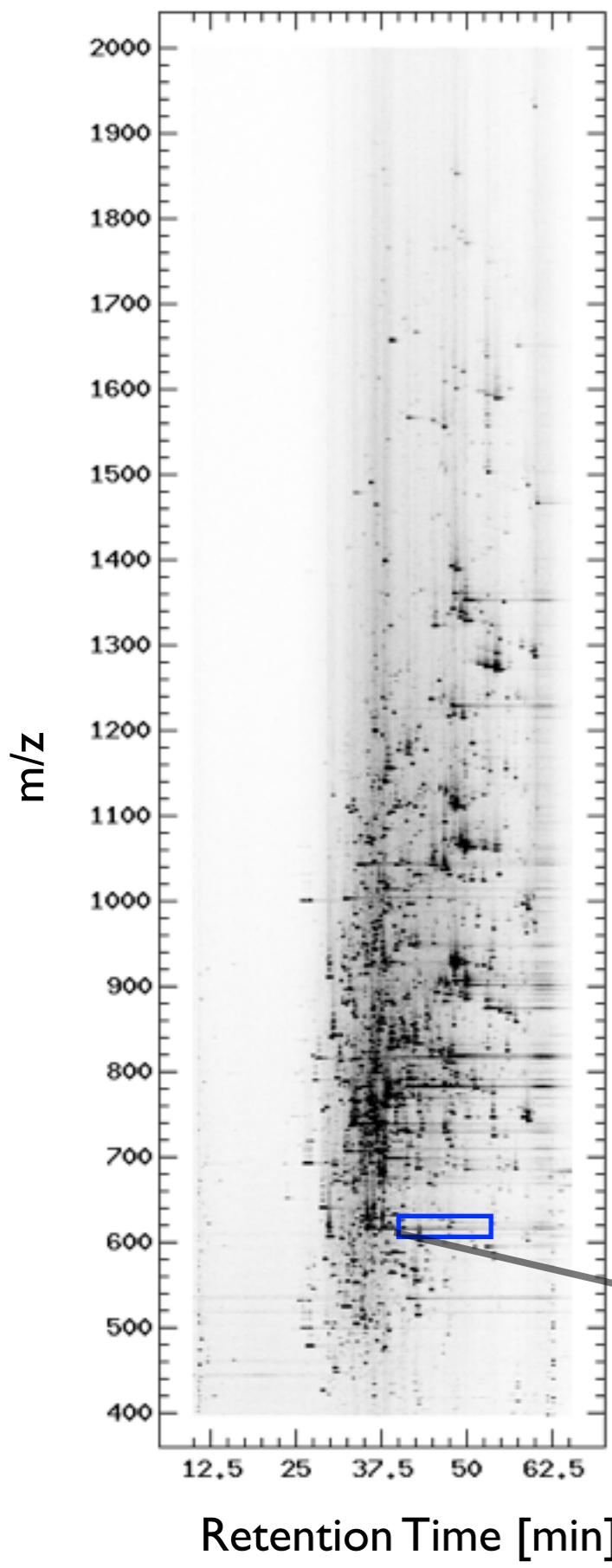


Extracted Ion Chromatogram

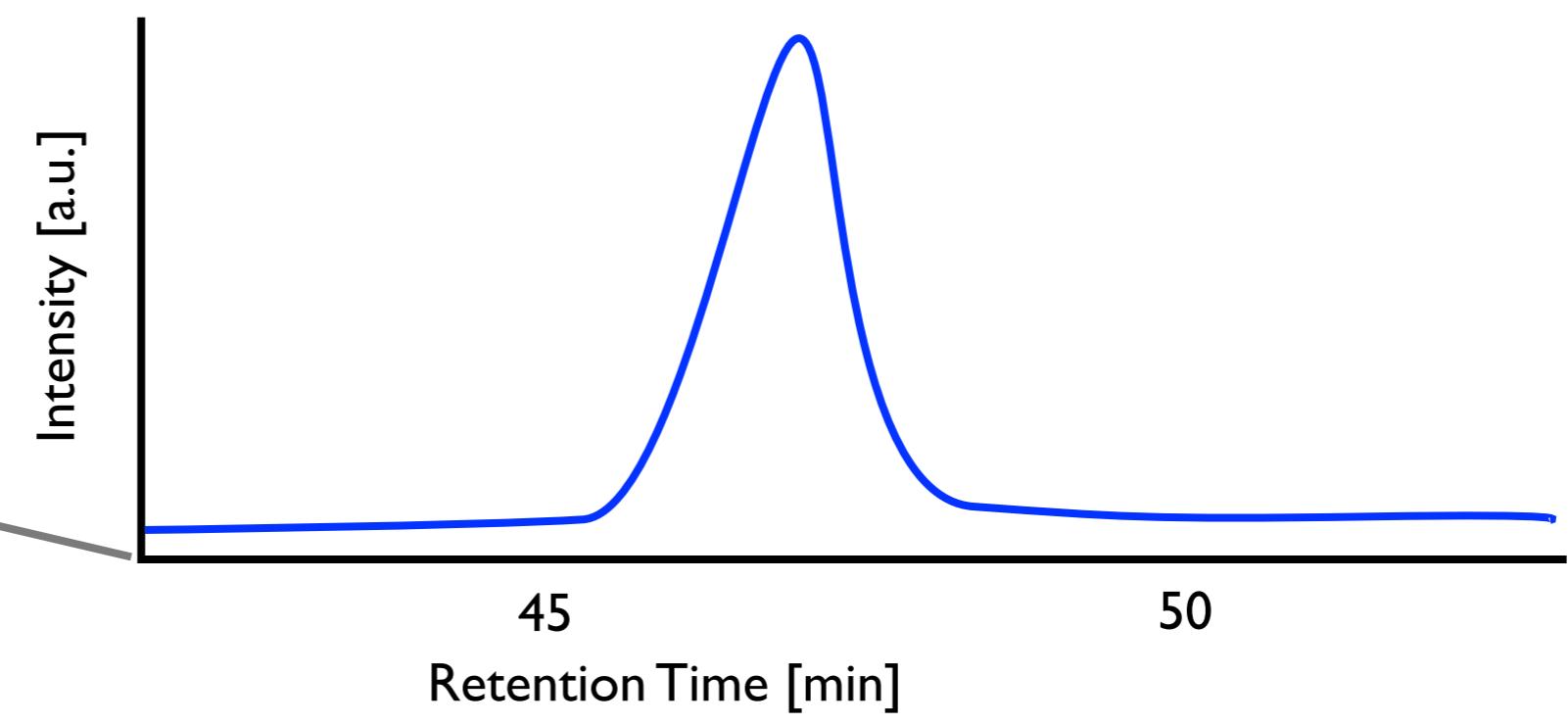


Chromatograms and Fragmentation spectra

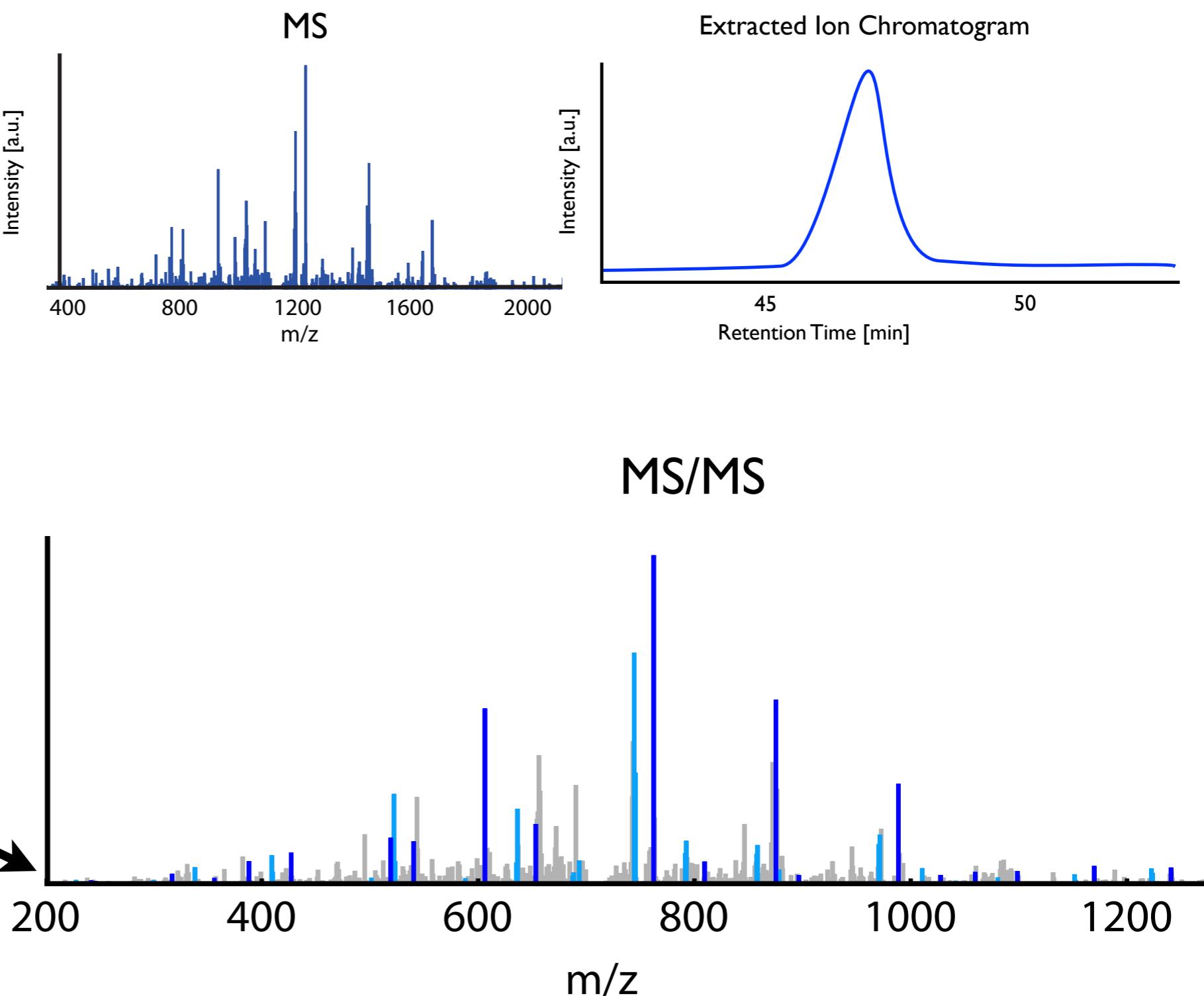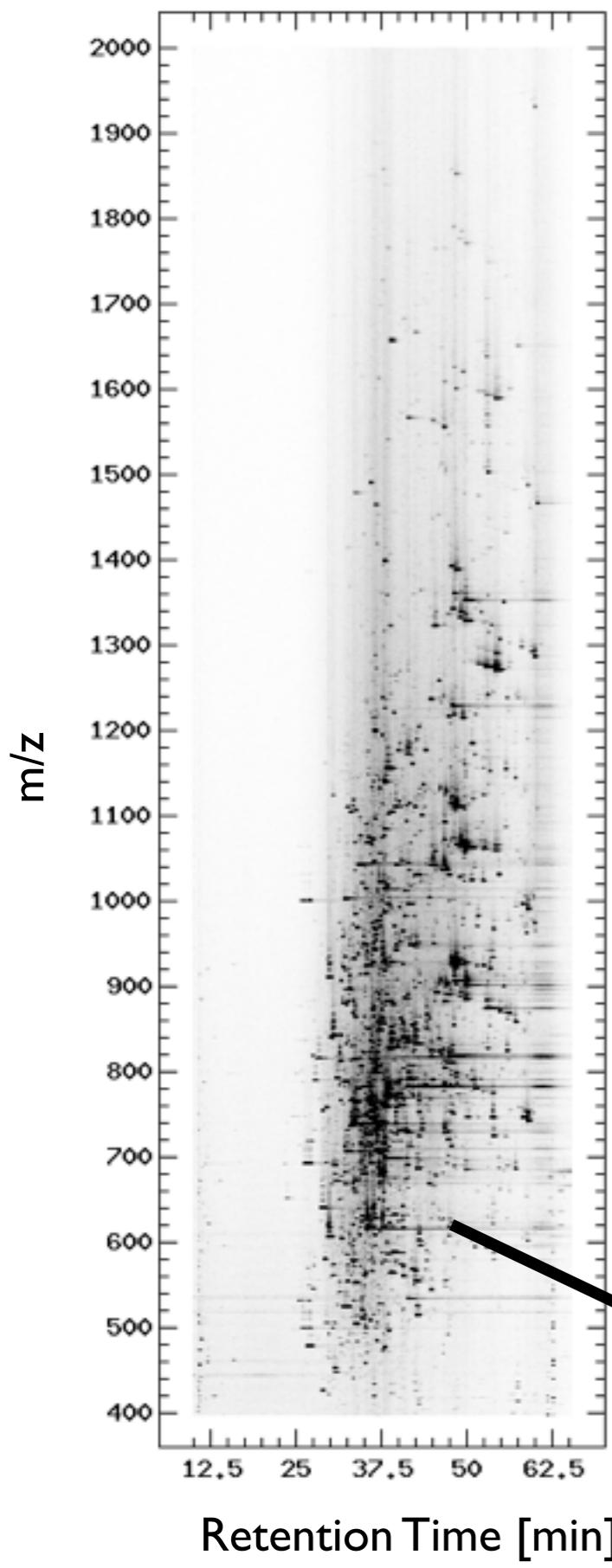

Peptide Identification

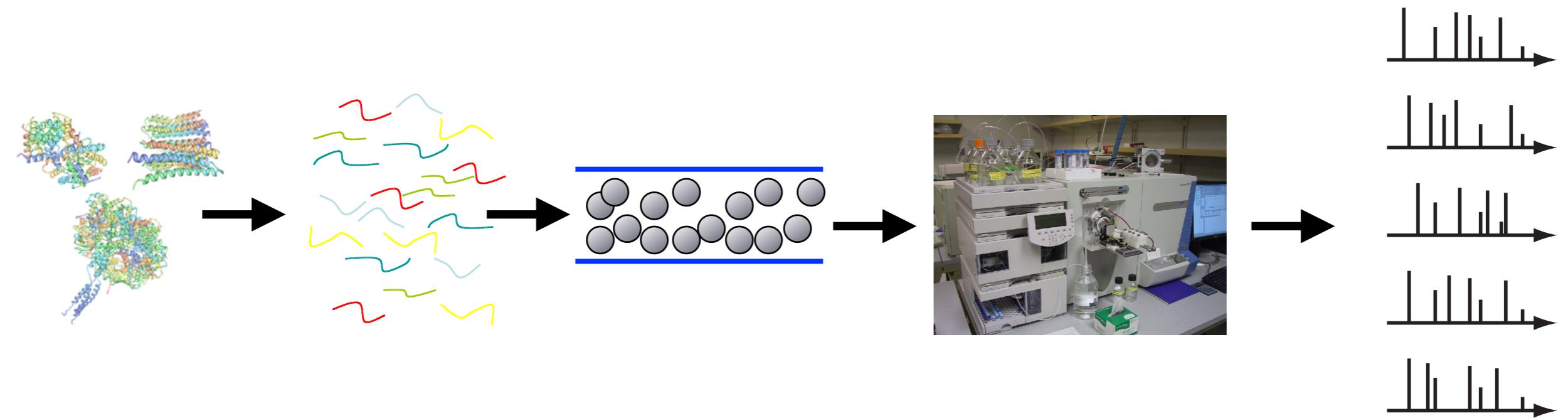


Peptide Identification

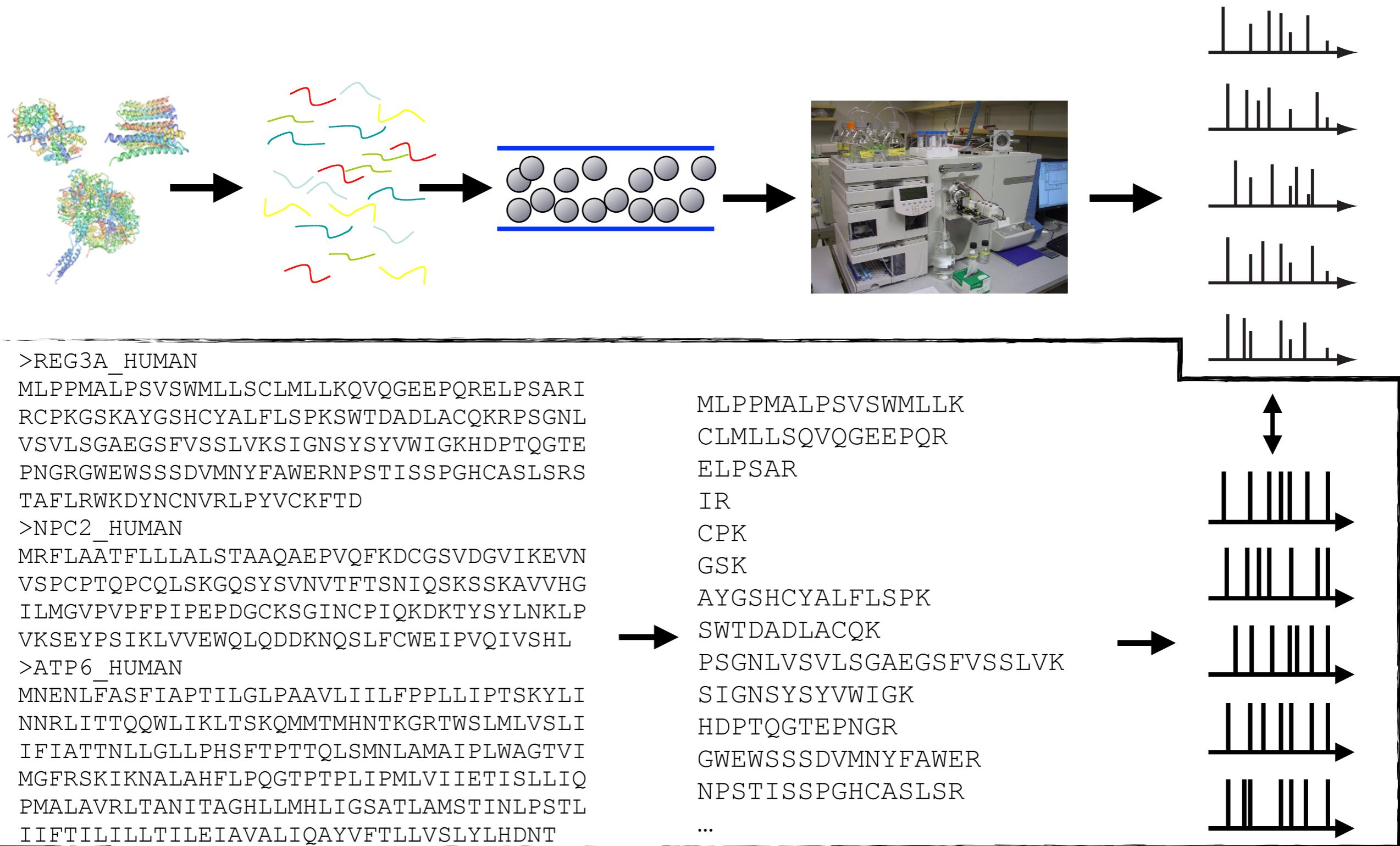

Quantitative Proteomics

Calibration curve

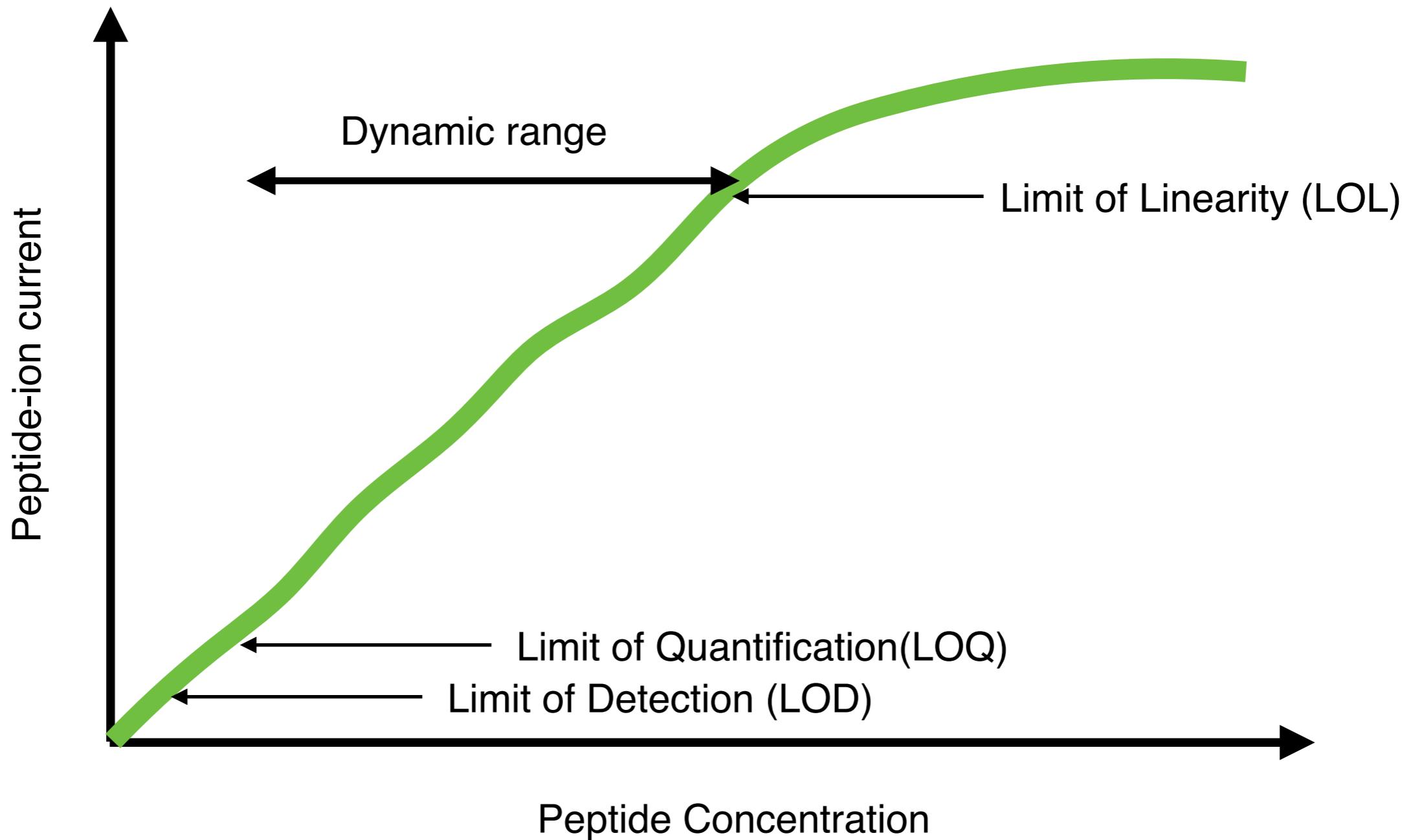


Different peptides have different slope

Why is protein quantification with mass spectrometry non-trivial?

Many of the steps from a protein to an ion current are non-predictively non-linear like:


- The efficiency of the enzyme
- The yield of the chromatography (some peptides do precipitate)
- The efficiency of the electrospray

Hence, relative quantification by mass spectrometry is easier than absolute quantification

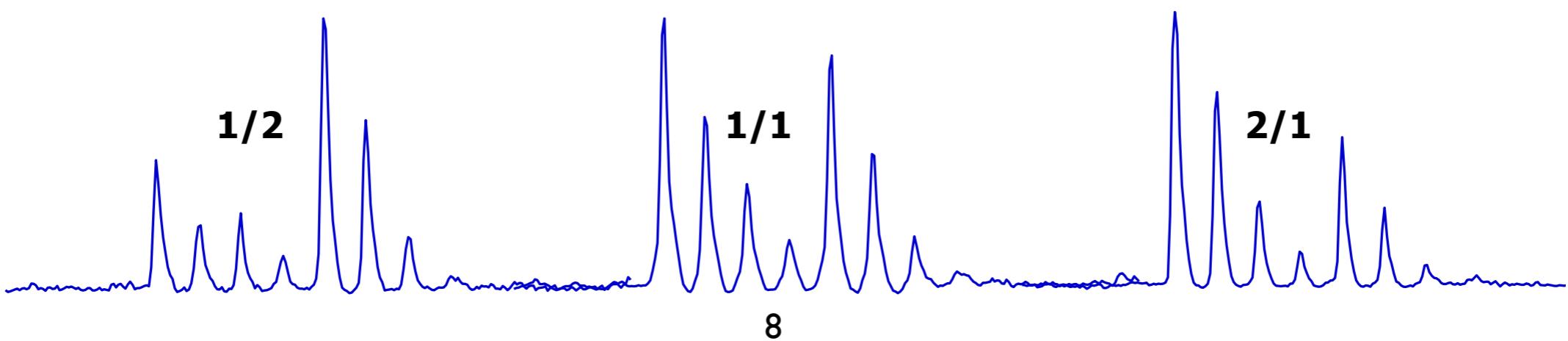

The primary principles in quantitation

- Make each sample distinguishable by either:
 - (I) introduce mass differences between the samples
 - (2) perform distinct experimental runs for each sample

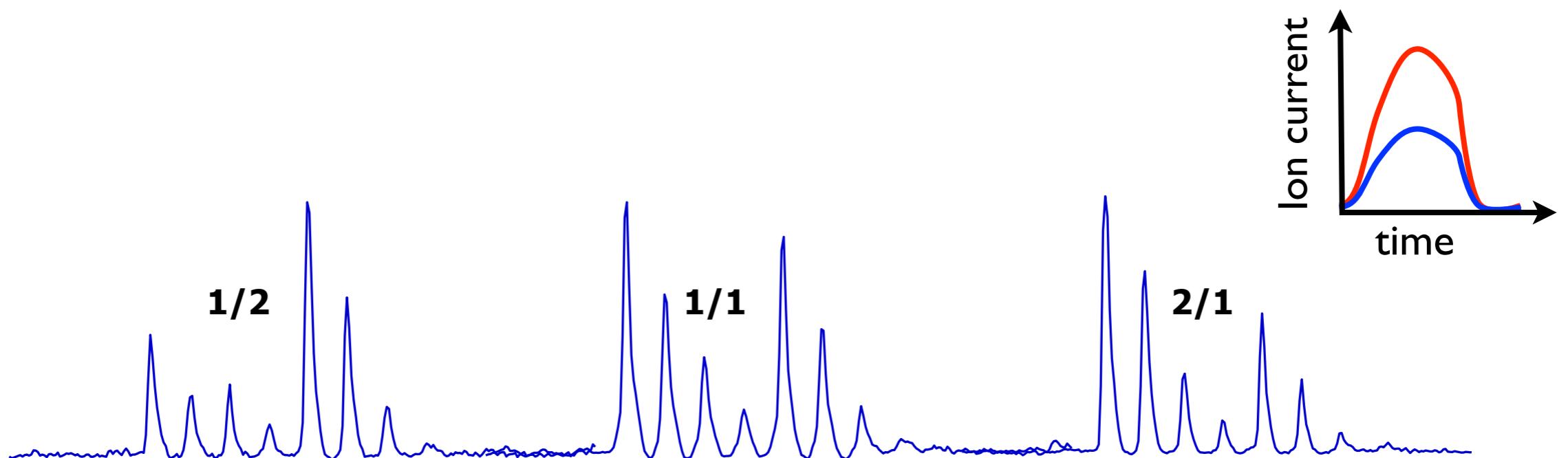

Protein quantitation

LC-MS(/MS) - protein or culture level labeling

- SILAC (2), cell cultures, relative

LC-MS/MS - peptide level labeling


- iTRAQ (2), peptides, relative
- Spiked peptides (eg. AQUA) (2), peptides, absolute

LC-MS(/MS) - label free

- MS Alignments (3), peptides, relative or absolute
- Spectral counting (3), peptide fragments, relative


Selected Reaction Monitoring

- With or without labels

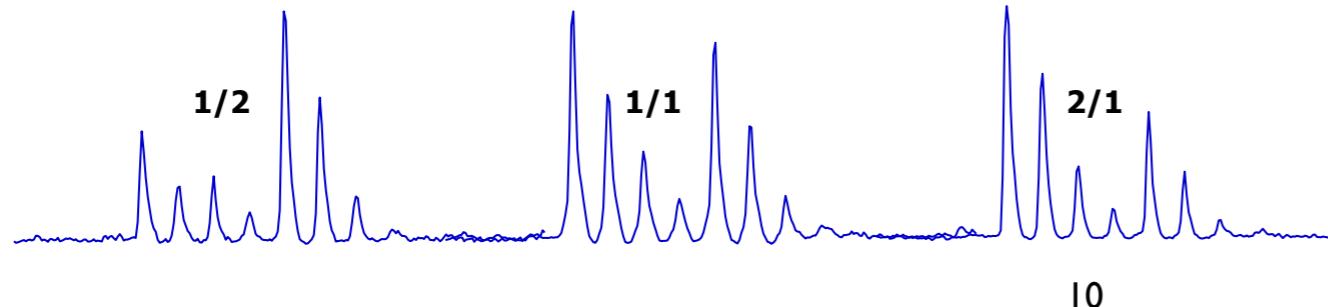


Culture-level labeling

SILAC (in-vivo heavy isotope label)

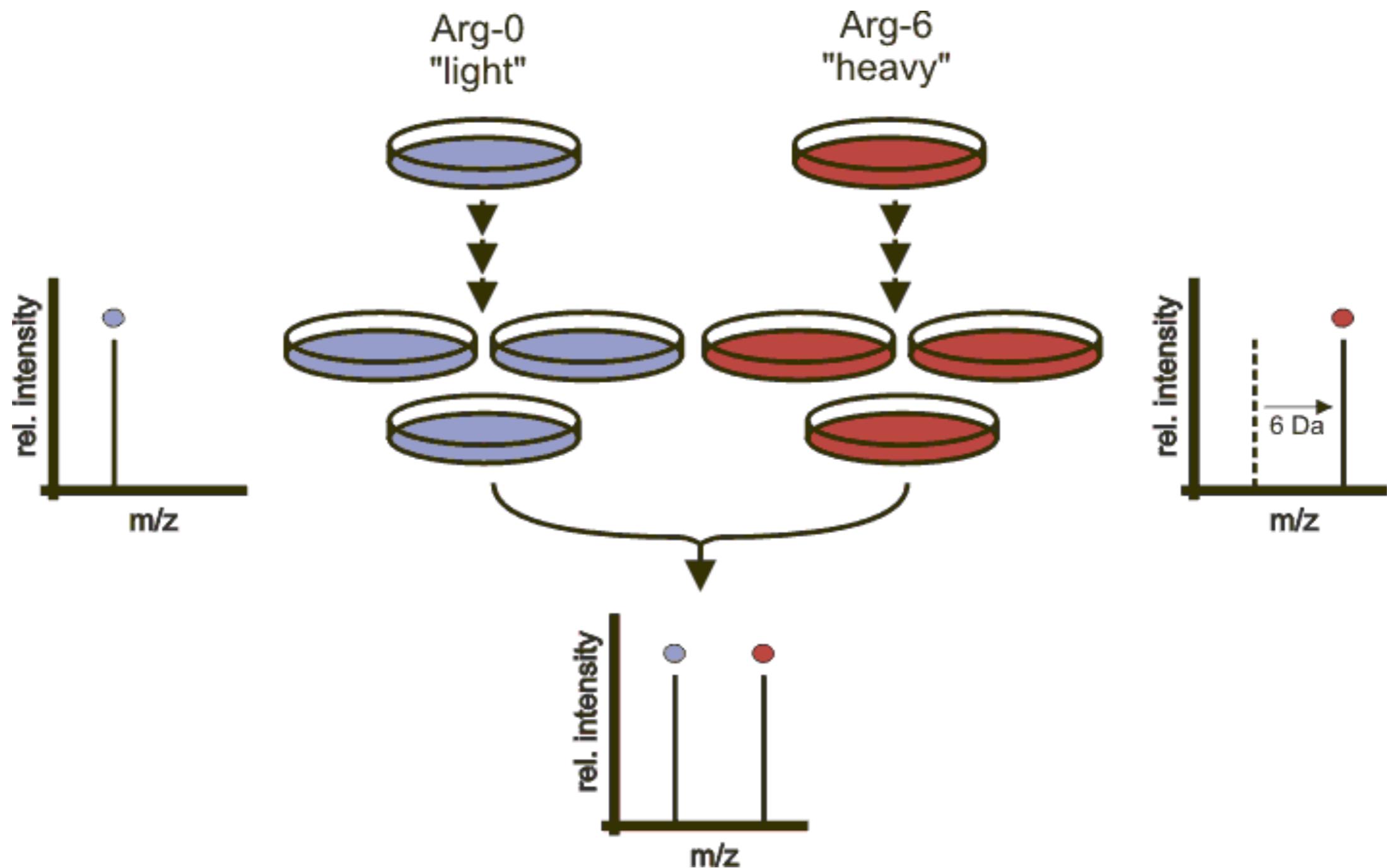



SILAC Mouse

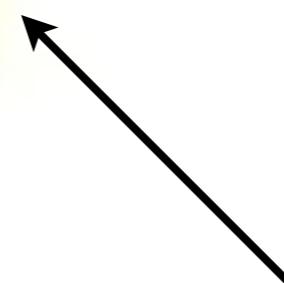

Heavy Cheese


Peptide-level chemical modifications

iTRAQ - MS/MS fragment (isobaric) labeling of peptides

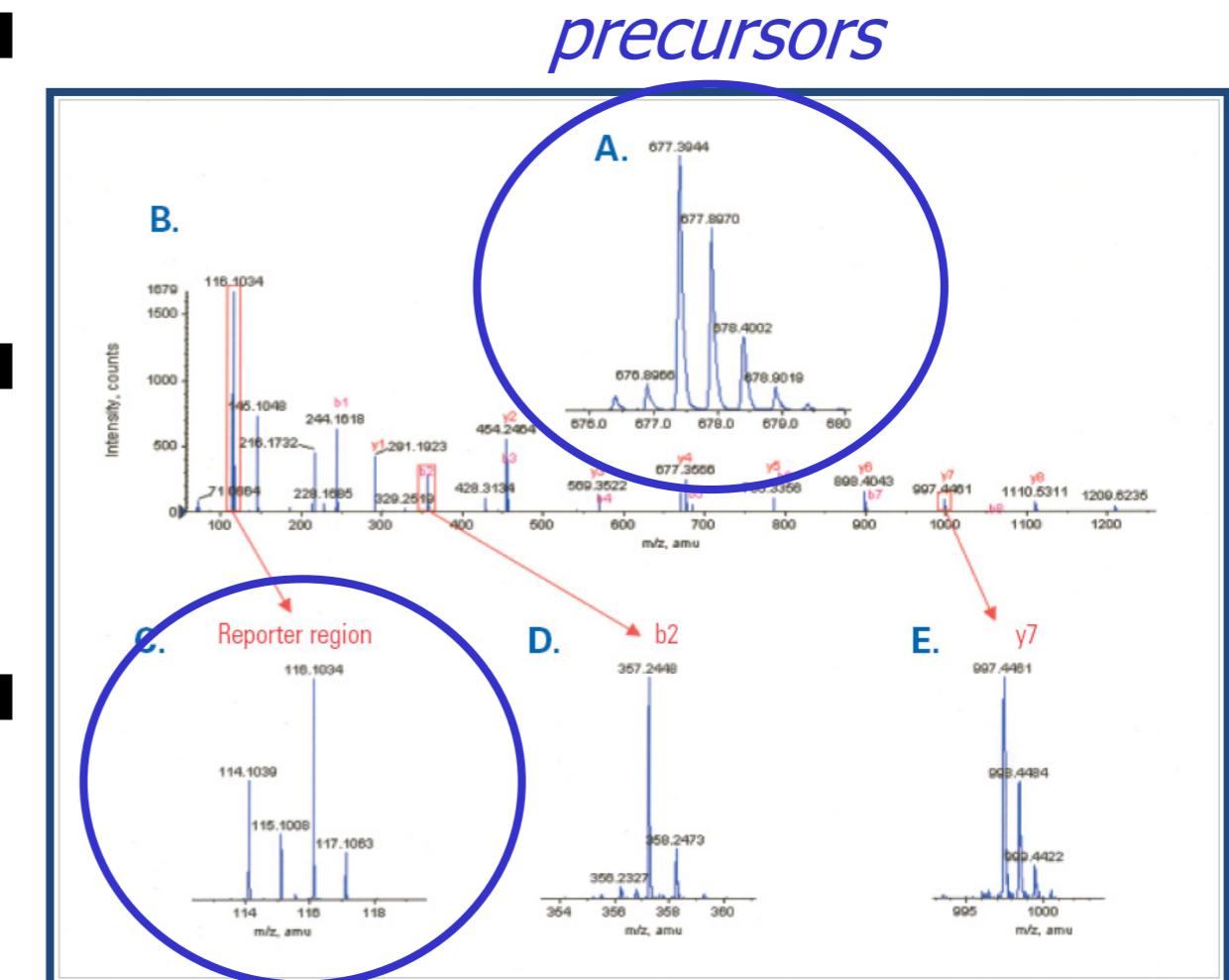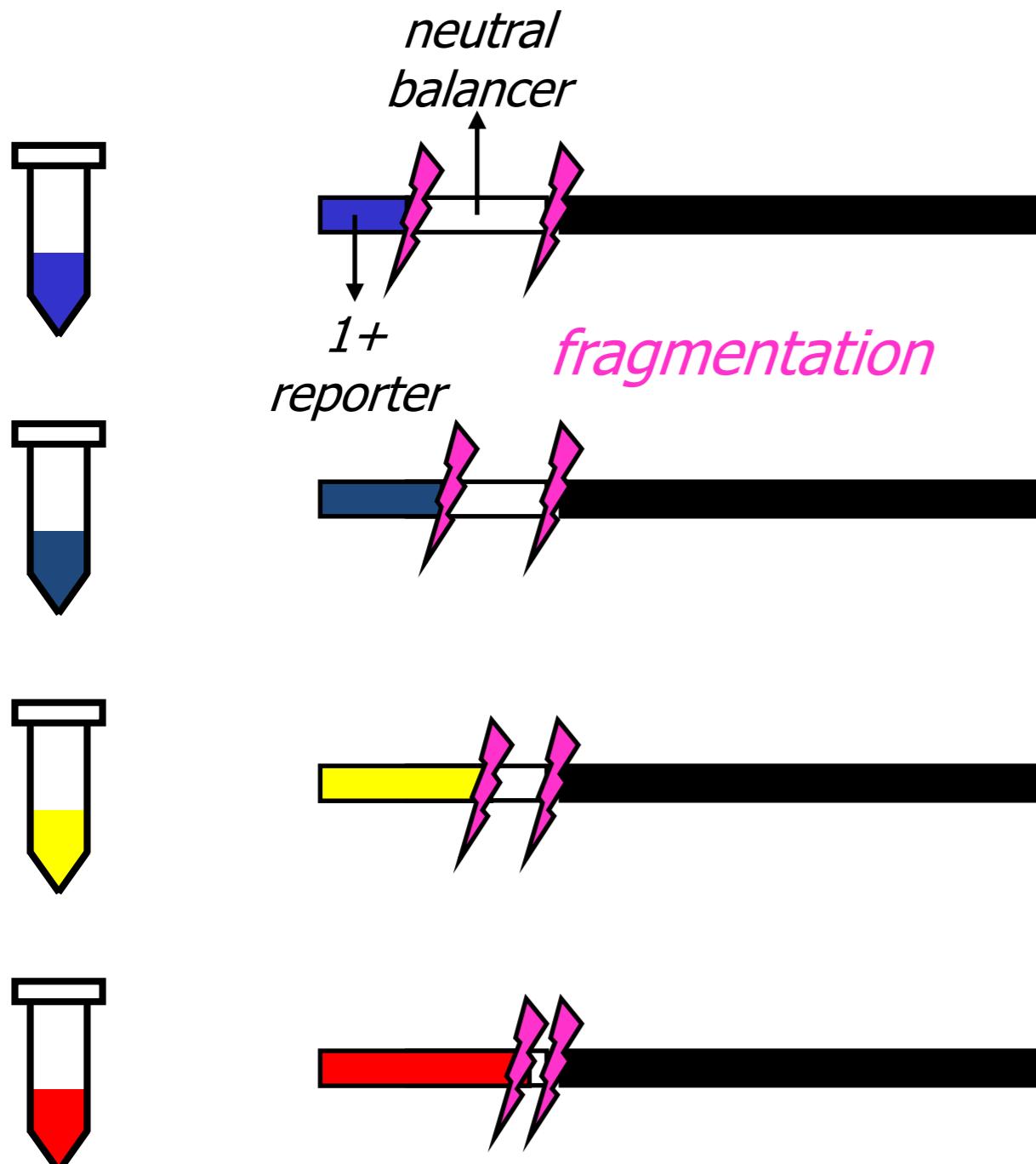


Relative sample abundance visible in fragmentation spectra

TMT - Yet another iso-baric label

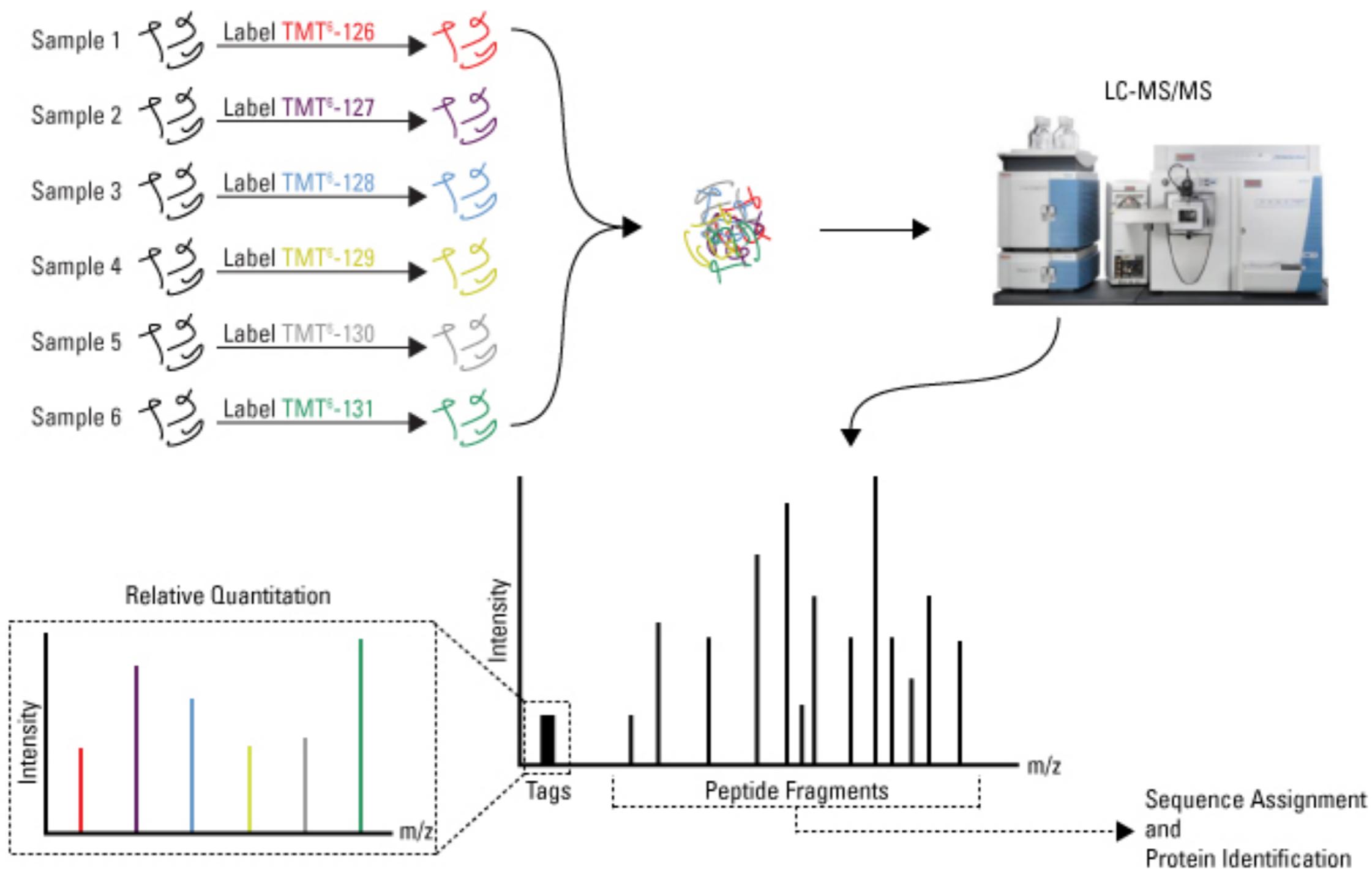


AQUA (peptide heavy isotope label)

Aimed at absolute quantitation

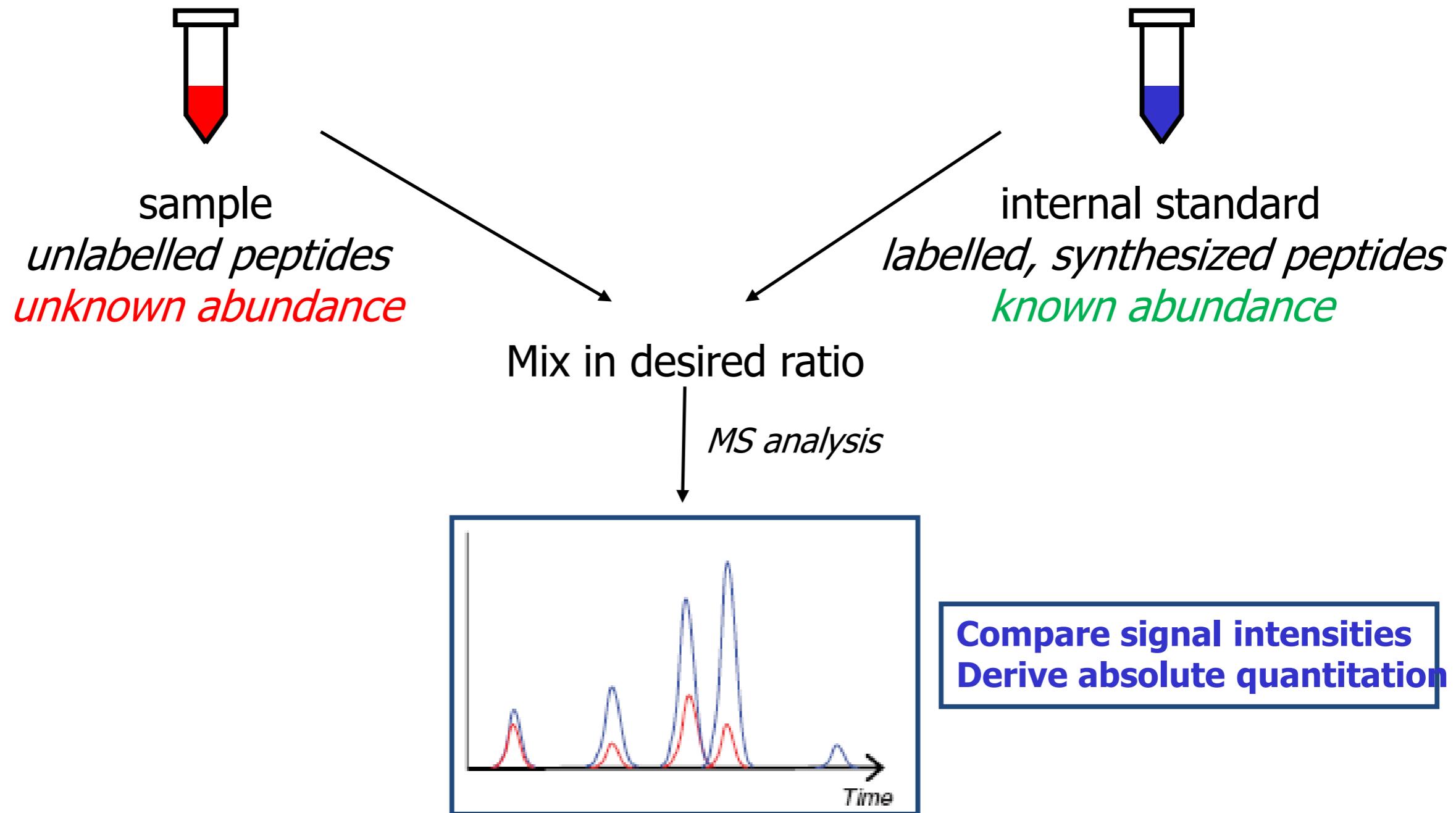


Labeled Quantification methods

At what point are the samples combined?

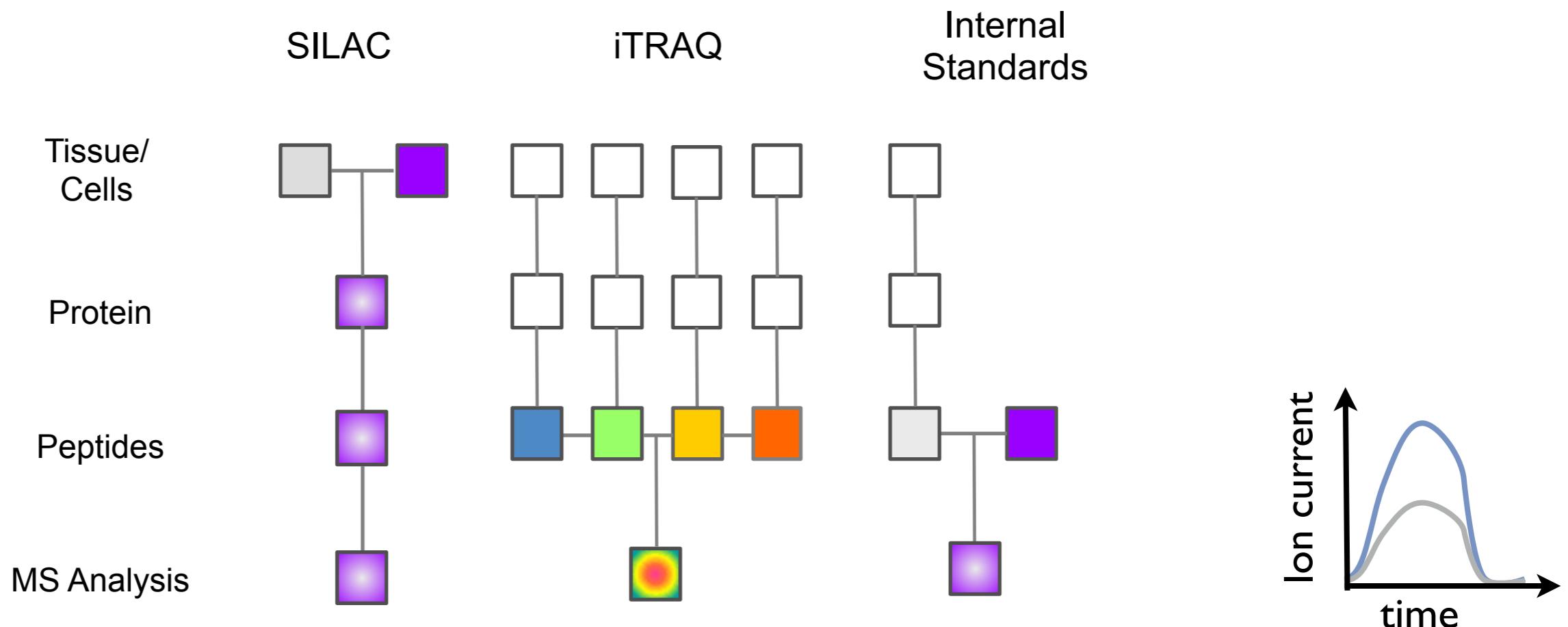

Label Free Quantification

Spectral counting

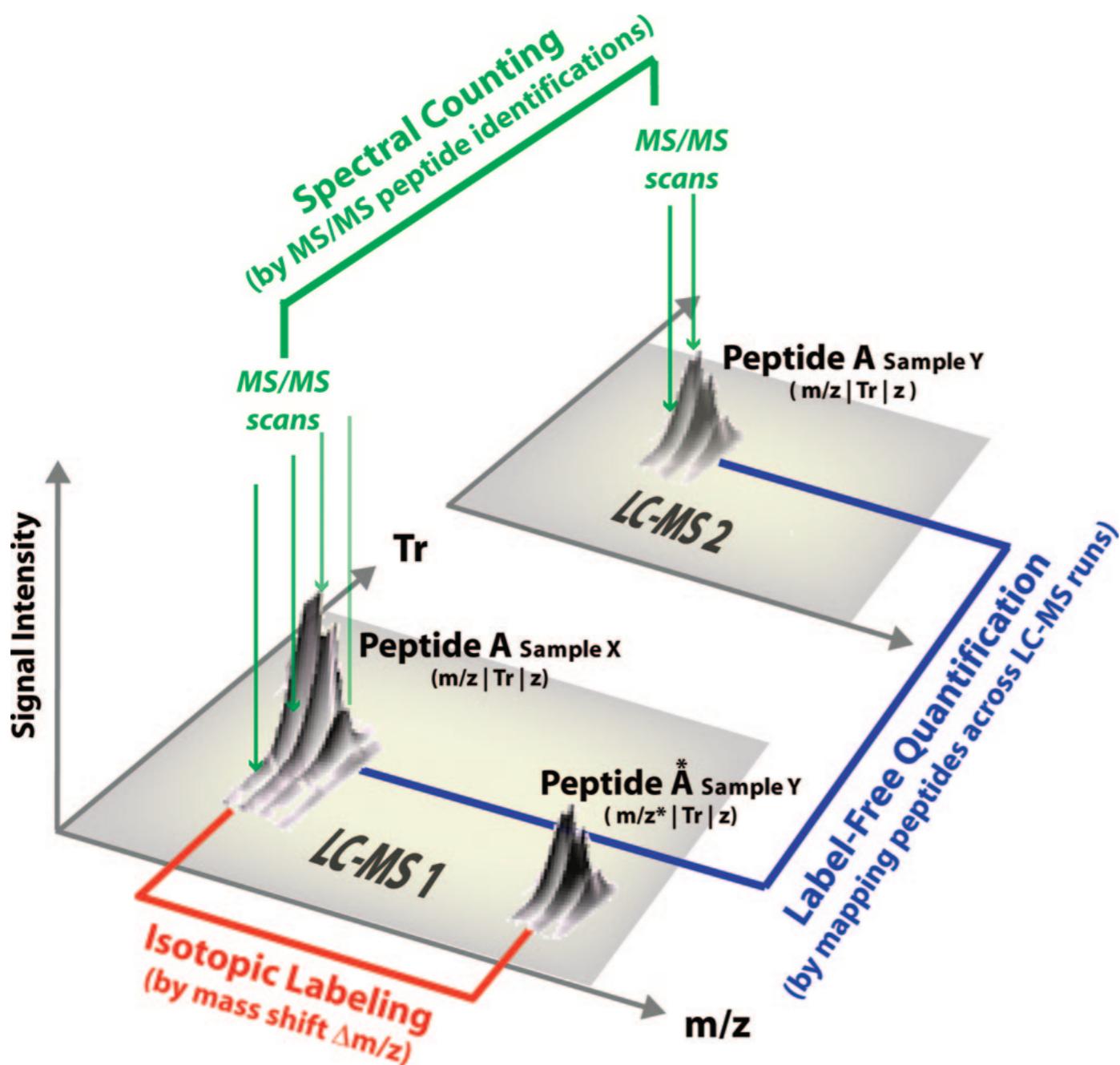

- (I) Count the number of fragmentation events that are triggered


emPAI (Exponentially Modified Protein Abundance Index)

$$\text{emPAI} = 10^{\text{PAI}} - 1$$

Where $\text{PAI} = \frac{N_{\text{observed}}}{N_{\text{observable}}}$



What is an 'observable' peptide?

- Peptides with a precursor mass between 800-2400Da.
- There is a roughly linear relationship between log protein concentration and the ratio of 'observable' peptides observed in range of 3-500 fmoles.
- If you know how much total protein you analyzed you can derive absolute abundances.


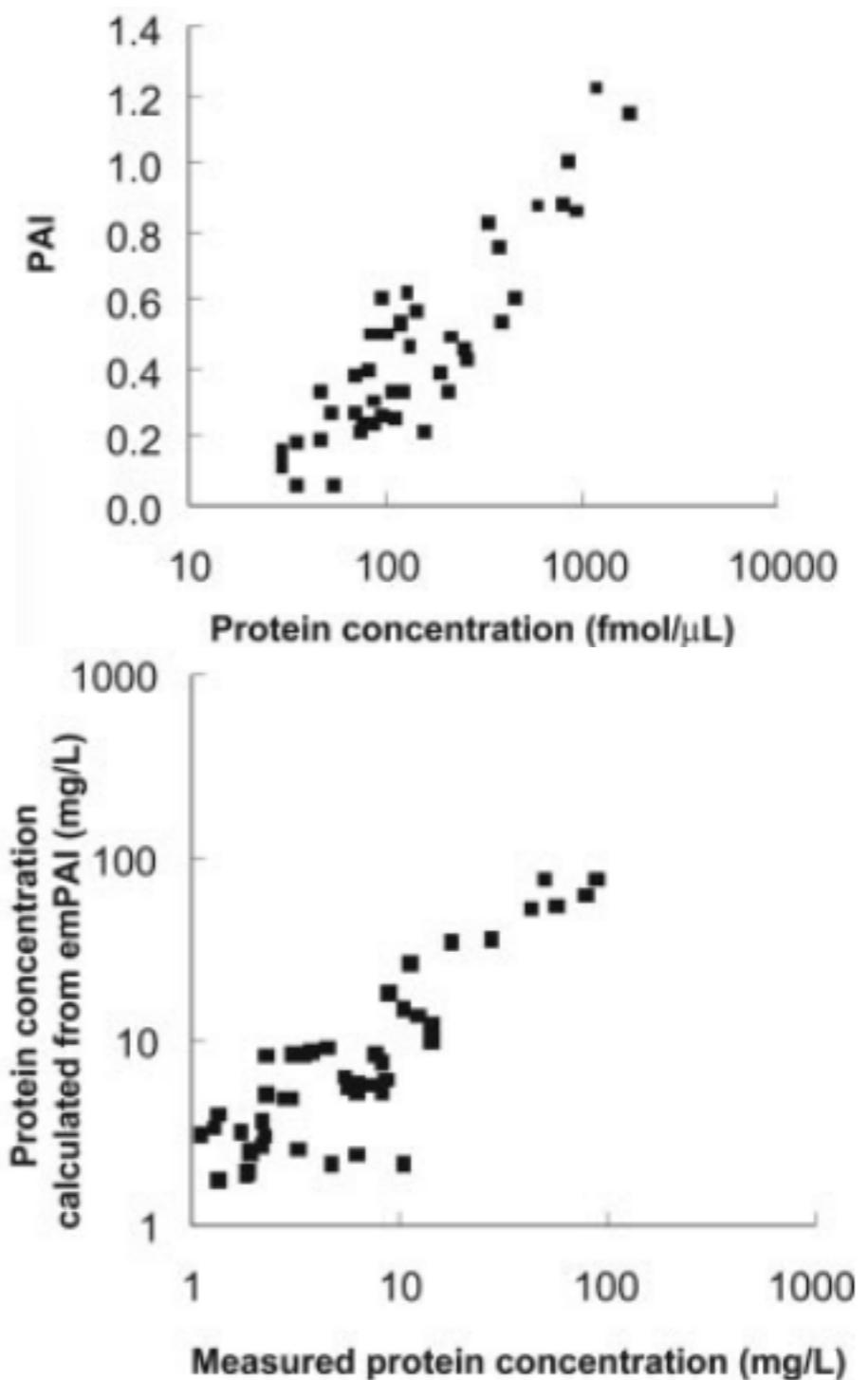


XIC (Extracted Ion Chromatogram)

- Measure intensity of peak during its elution off HPLC column and into the mass spectrometer.
- Measure area of peak in XIC.
- More accurate than selecting peak intensity for one given scan.

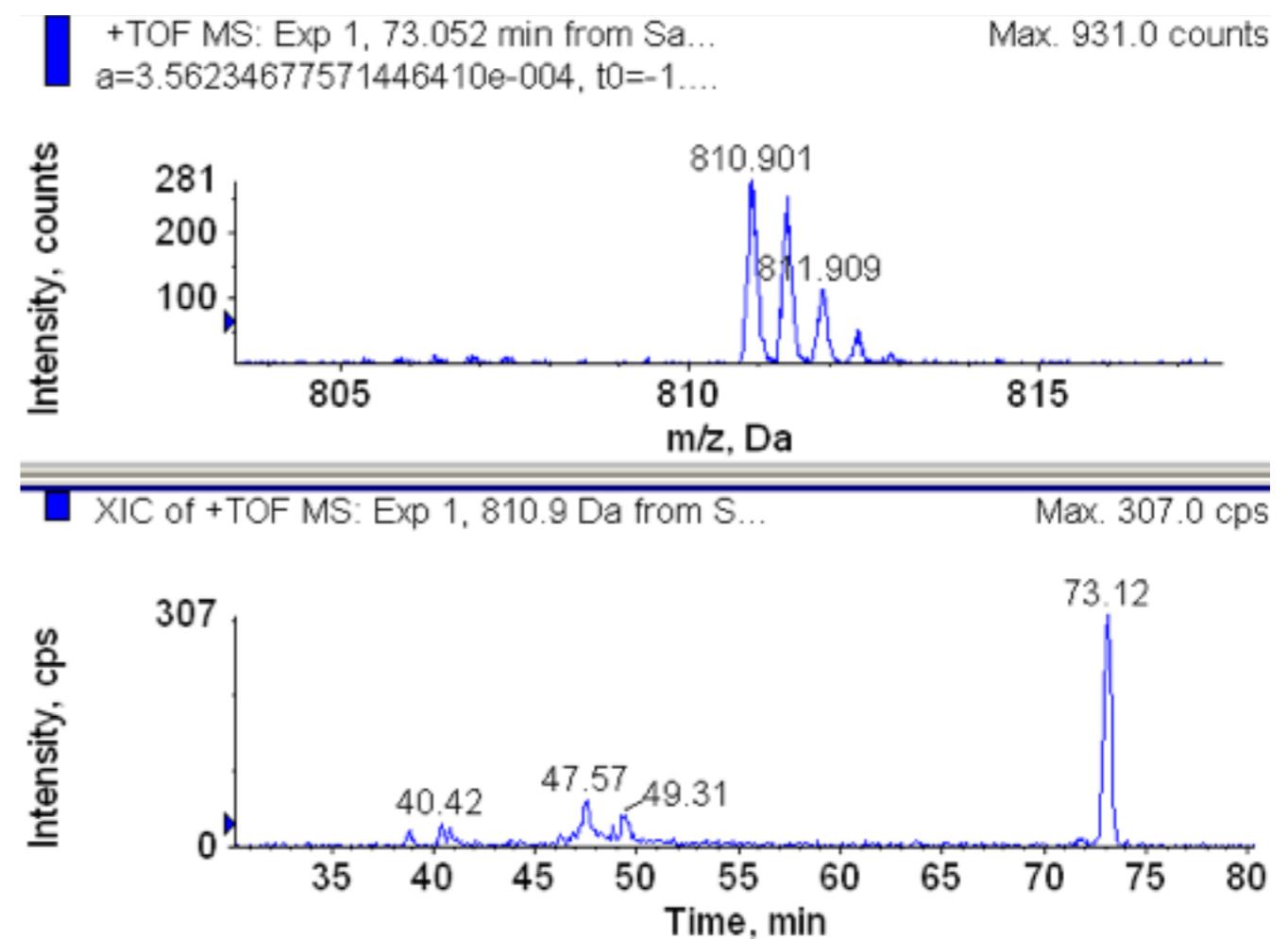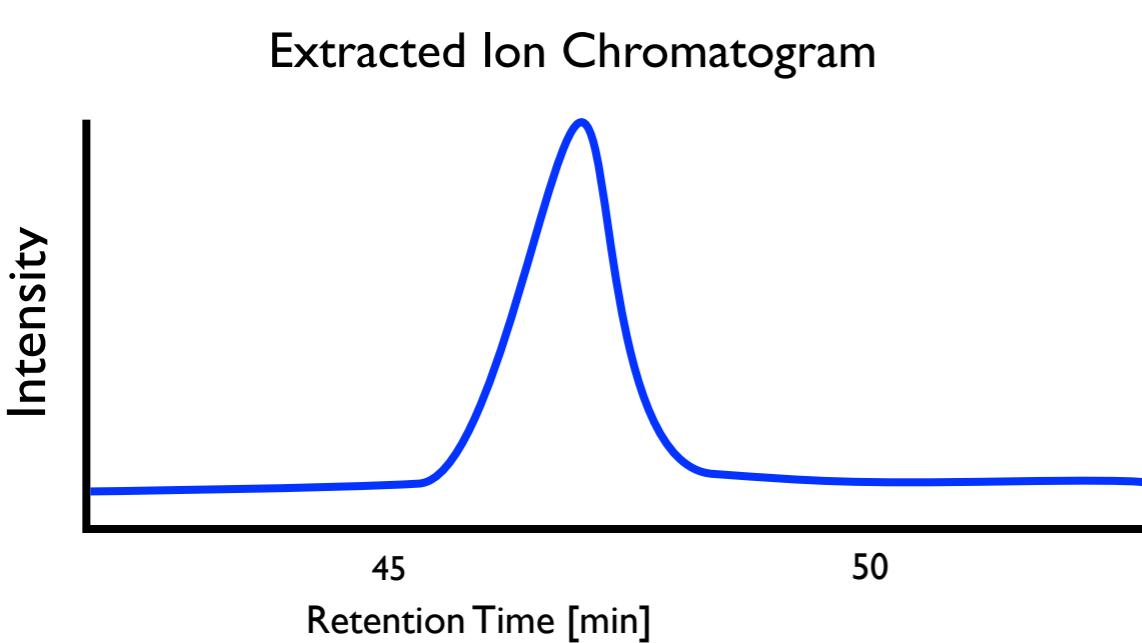


Label Free: MS Alignments

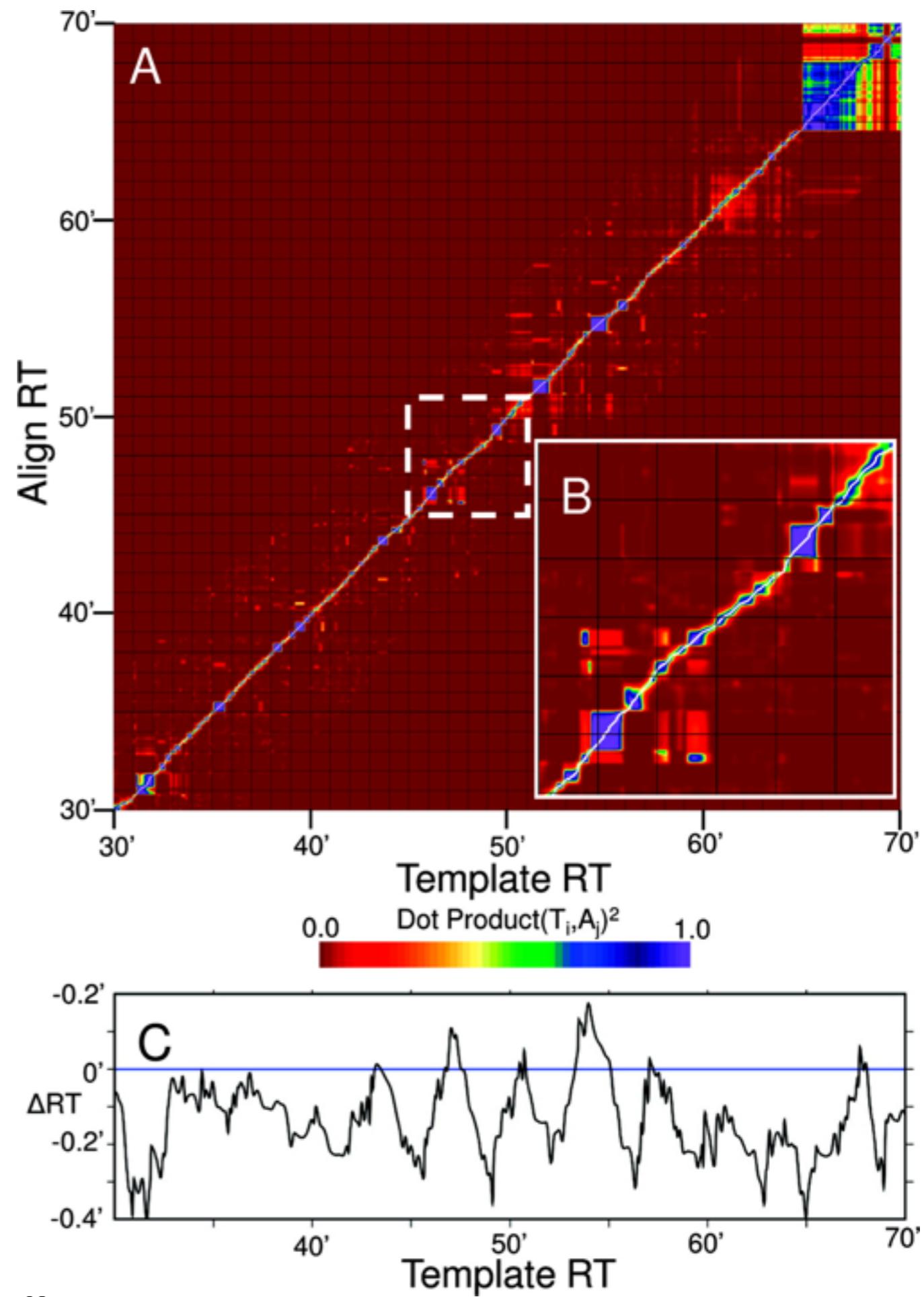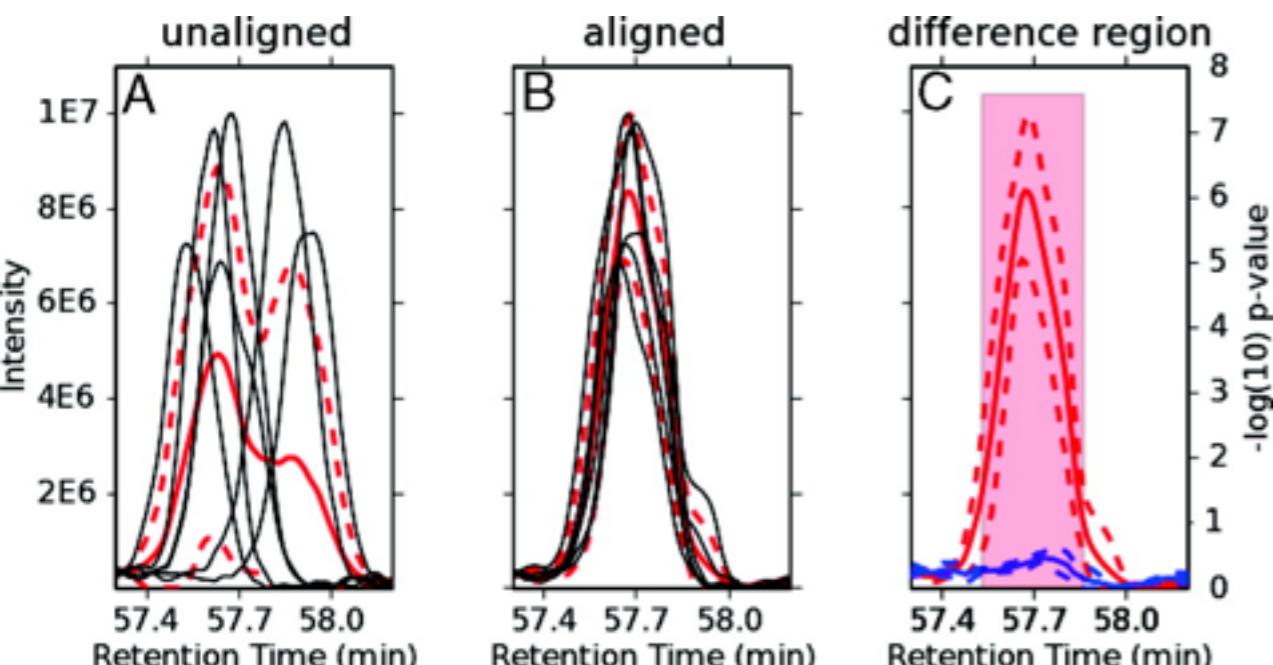

[Finney et al JPR 2008]


Other Acquisition Techniques

Selected Reaction Monitoring (SRM or MRM)

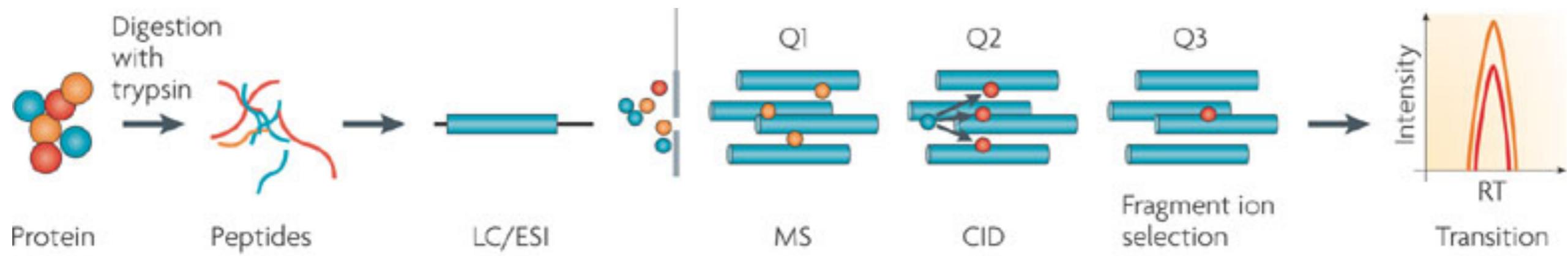

- + Higher sensitivity (low-attomolar)
- + Broader dynamic range (up to 5 orders of magnitude)
- Requires known protein targets and their characteristics
- max 200-300 measurements per run

Data Independent Analysis (DIA or SWATH)

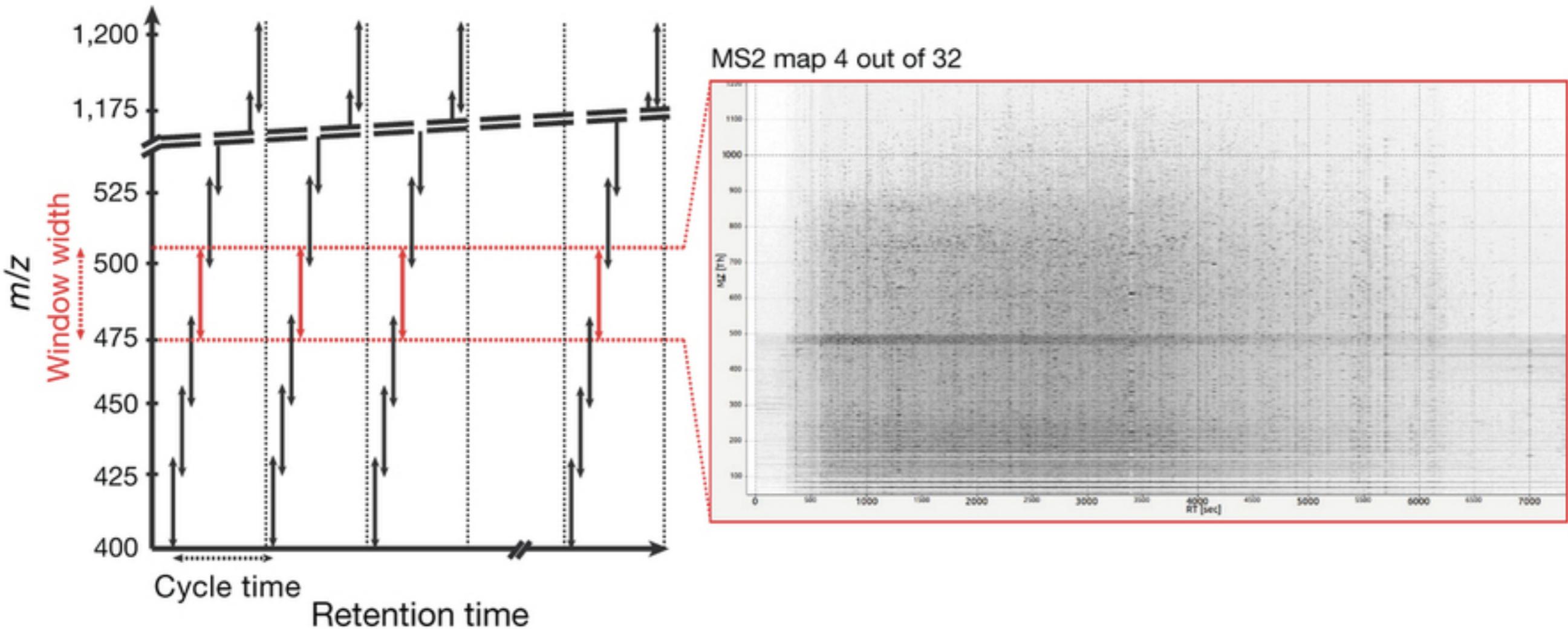

- + Highly reproducible
- + Broad dynamic range (~4 orders of magnitude)
- + Measures a large number of proteins per run
- Requires known protein targets and their characteristics

Comparison of acquisition techniques


	Shotgun	SRM	DIA
Quantifiable Proteins per run	1000s	100s	1000s
Reproducibility	++	++++	+++
Sensitivity	+++	++++	+++
Dynamic Quantification Range	3-4 magnitudes	5 magnitudes	3-4 magnitudes

Selected Applications

Identification of components in protein complexes


Tandem Affinity Purification

- An 'affinity tag' is first attached to a target protein (the 'bait')
- Bait proteins are precipitated, along with its protein complex
- Purified protein complexes are resolved by one-dimensional SDS-PAGE
- Proteins are excised from the gel, digested with trypsin, and analyzed by mass spectrometry.

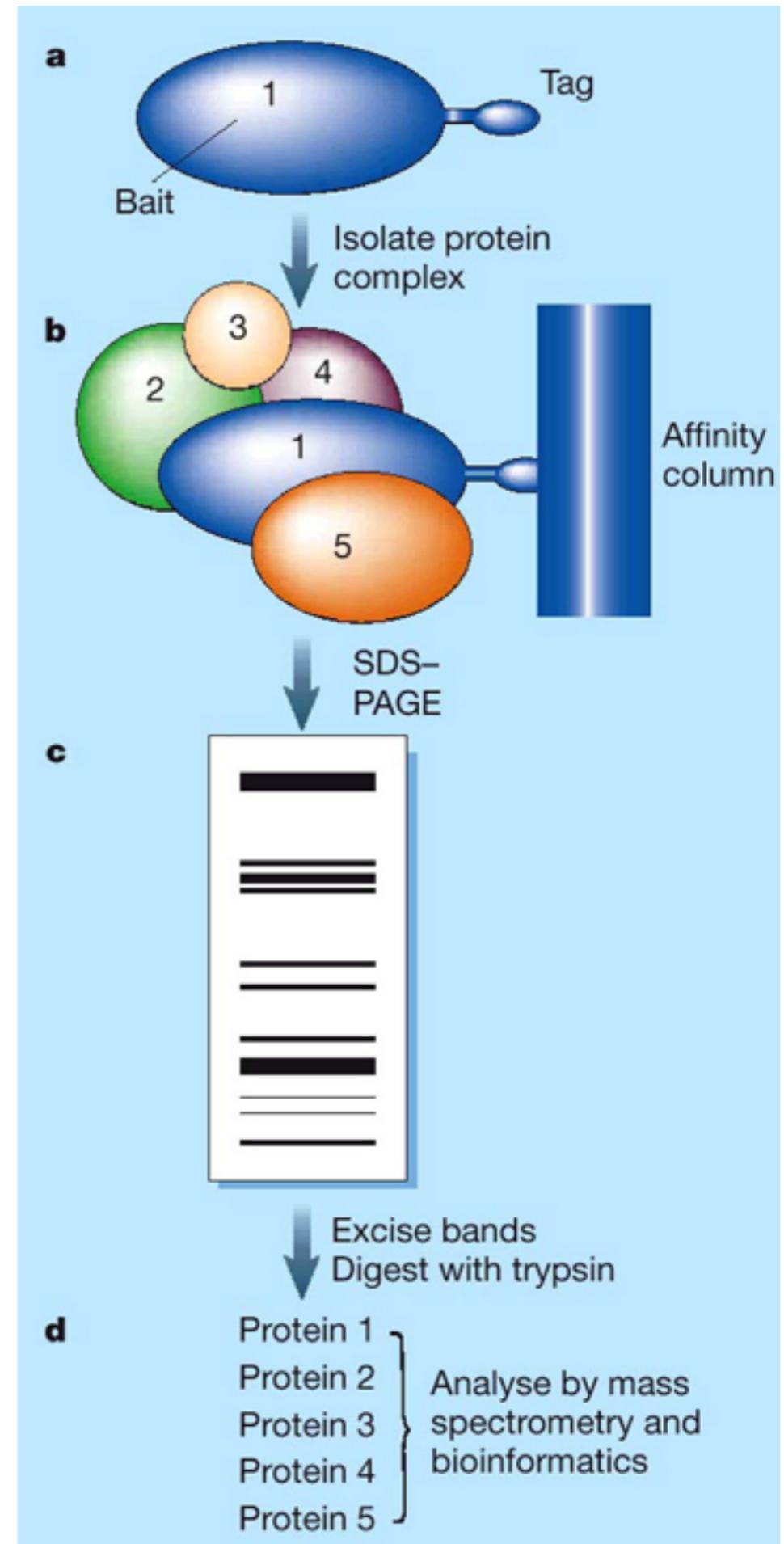

[Kumar & Snyder, 2002]

Example: Barcoding cell types with Cell surface–capturing (CSC)

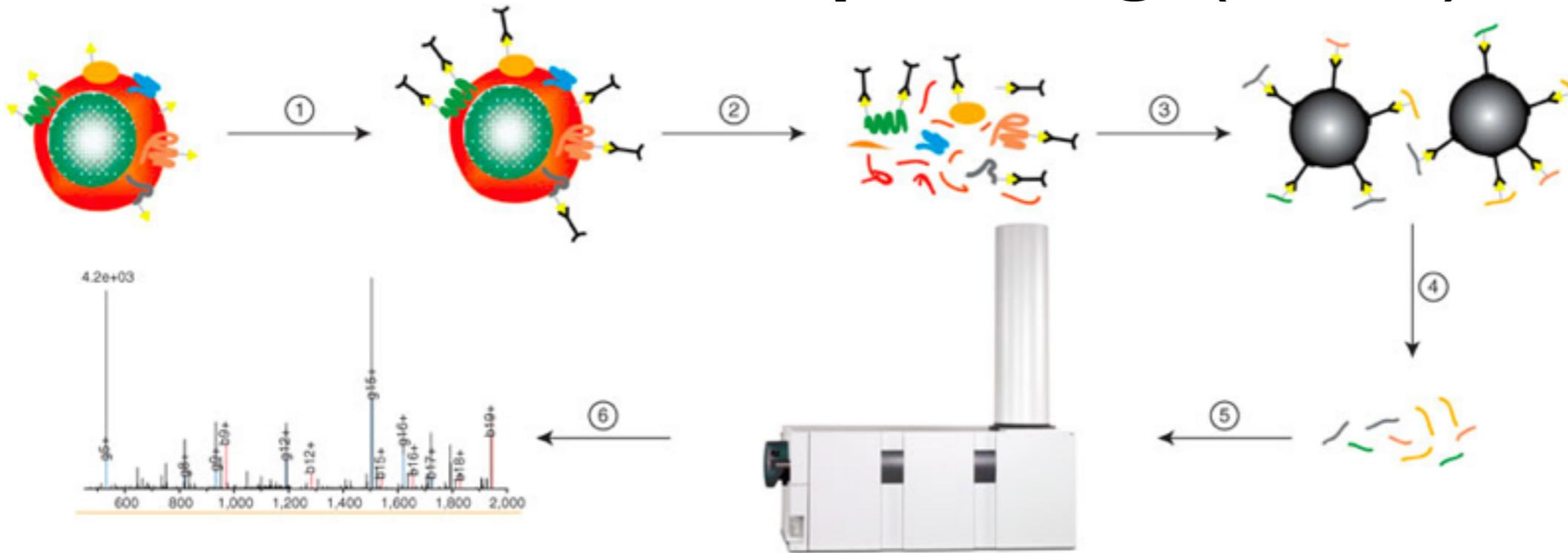

Identifies 110 cell surface proteins

Proof of principle: differentiate
Ramos B cells from Jurkat T cells

[Wollscheid et al., Nat Biotech 2009]

